

ZOLL AED 3

Accessories and Training Products

ZOLL®



ACCESSORIES TO SUPPORT YOUR RESCUE


CONSUMABLE RESCUE ACCESSORIES

Long-lasting for Every Rescue

CPR Uni-padz® III electrodes

Truly universal, these pads can be used in an adult or child rescue. For child rescue, reposition the CPR Uni-padz III as shown and activate the child rescue mode. CPR Uni-padz III electrodes provide full CPR feedback for adult patients and have a five-year installed shelf life.

Part # 8900-000280-01

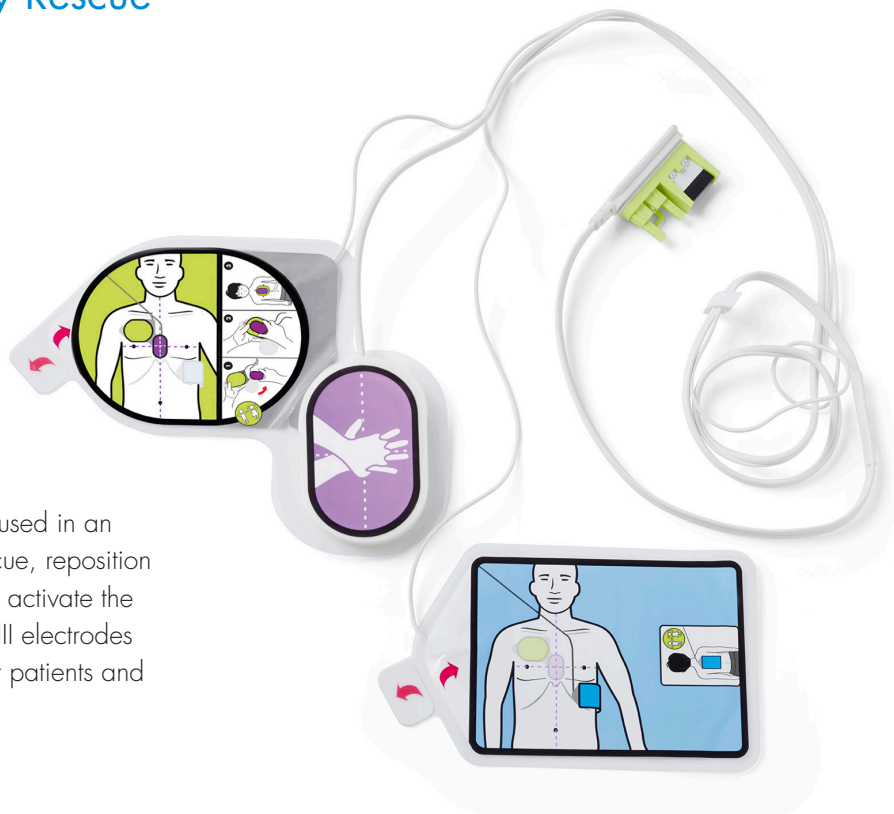

Every set of CPR Uni-padz electrodes comes with a Rescue Accessory Kit that includes scissors for cutting clothing to expose the victim's chest, non-latex gloves and a breathing mask for rescuer protection, paper towel to ensure dry contact with electrodes, razor for removing hair, and a wipe for cleanup.

Smart Battery

The ZOLL AED 3 uses a smart battery that can report its power level status on the device's LCD screen or transmit it automatically over WiFi to your AED program management solution. In stand-by mode, the ZOLL AED 3 disposable lithium manganese dioxide battery has an installed life of up to five years.

Part # 8000-000696



ZOLL AED 3 CLINICAL UNIT TRAINING ACCESSORIES

Tools to Simulate Every Scenario



Training/Demonstration CPR Uni-padz Electrodes

Connects to a ZOLL AED 3® defibrillator for training and simulations. Includes one set of two reusable training pads
Part # 8900-000284



ZOLL AED Simulator

Generate either shockable or non-shockable heart rhythms for various rescue scenarios. Simulator uses a 9-volt battery, included with purchase.
Part # 8000-000925



CPR Uni-padz Training Replacement Pads

Eight sets of two reusable training pads for use for simulation and training
Part # 8900-000285



Demo Manikin

Adult manikin for AED/CPR training simulates a victim of 130 pounds (60 kg). Removable head for compact storage and hands-only training
Part # 8000-0835-01

CPR Uni-padz Travel Trainer

One-piece, compact, easy-to-carry demonstration tool simulates shockable and non-shockable heart rhythms.
Part # 8008-000053-01





AED SIGNAGE

Create High Visibility

AED Sign Package

3D wall sign and 8.5" x 11" wall-mounted sign to improve AED visibility

Part# 8000-0825



Strobe for Standard Surface-mounted Wall Cabinet

Part# 8000-001259



Strobe for Recessed Wall Cabinet

Part# 8000-001267



WALL CABINETRY

Secure & Accessible

Keep your ZOLL AED 3 defibrillator secure and easily accessible for medical emergencies with a variety of battery-operated, alarmed wall cabinets.* Each holds a ZOLL AED 3, a backup set of electrodes, and a backup

battery should you need them. Choose the surface-mounted model for easy installation or the semirecessed or fully recessed models for less intrusive storage. All three models include installation instructions.



Semirecessed Wall Cabinet

Part# 8000-001257

Meets ADA (Americans with Disabilities Act) requirements
14 inches x 14 inches x 3.5 inches
(35.5 cm x 35.5 cm x 8.9 cm)



Fully Recessed Wall Cabinet

Part# 8000-001258

Meets ADA (Americans with Disabilities Act) requirements
14 inches x 14 inches x 1.5 inches
(35.5 cm x 35.5 cm x 3.8 cm)



Standard Surface-mounted Wall Cabinet

Part# 8000-001256

17.5 inches x 17.5 inches x 7 inches
(45.5 cm x 45.5 cm x 17.8 cm)

*9-volt battery included

WALL MOUNTS

Simple & Universal

The universal wall mount holds the ZOLL AED 3 alone or in its carry case. It includes mounting hardware, consumable tracking card, wall-mount template, and installation instructions.



Wall Mount Bracket for ZOLL AED 3 (without carry bag)

Part# 8000-001255



Wall Mount Bracket for ZOLL AED 3 with carry case and backup CPR Uni-padz electrodes (in carry case)

Part# 8000-001266

CARRY CASES

Rugged & Portable

Choose from three different carry case models to safely carry the ZOLL AED 3 defibrillator: the standard carry case, the small hard-shell carry case, or the large hard-shell carry case. The hard-shell cases are waterproof, with an ingress protection (IP) rating of IP68 for dust and water.

Standard Carry Case

In addition to the ZOLL AED 3 with pre-connected CPR Uni-padz electrodes, the standard carry case can also hold a backup set of CPR Uni-padz.

Part# 8000-001250

Replacement Shoulder Strap For ZOLL AED 3 Carry Case

Part# 8000-001252



Large Hard-shell Carry Case

Protect and carry everything needed for a vigorous rescue with the large hard-shell carry case, including a backup set of CPR Uni-padz electrodes and a backup battery.

Part# 8000-001254



Small Hard-shell Carry Case

For lightweight maximum protection, choose the small hard-shell carry case. It accommodates the ZOLL AED 3 and one spare battery.

Part# 8000-001253



Battery Case

Store the ZOLL AED 3 battery in this convenient pouch for easy identification and access when needed. Battery case stores easily in both the large and small hard-shell carry case.


Part# 8000-001251

ZOLL AED 3 TRAINER AND ACCESSORIES

Simulate a Rescue with Real CPR Help

The ZOLL AED 3® Trainer is the only automated external defibrillator (AED) training solution that uses Real CPR Help® to provide real-time feedback on a trainee's compression quality.

- Trainees learn the proper rate and depth of CPR chest compressions.
 - An active LCD display with full-color rescue images mirrors what a student will see on a clinical ZOLL AED 3 device.
 - An on-device pause button allows the instructor to pause a training scenario.
 - Trainer device ships with training harness, training electrode gels/liner, batteries, and Operator's Guide.
- Part # 8028-000001-01

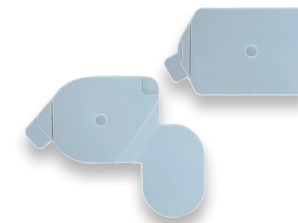

Replacement CPR Uni-padz Electrode Liners


Part # 8028-000013



CPR Uni-padz Electrode Training Harness

For use with ZOLL AED 3 Trainer

Part # 8028-000010

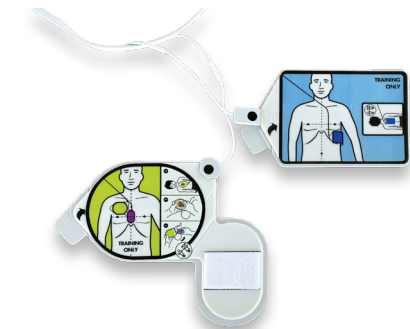

Replacement Adhesive Gels for CPR Uni-padz Electrodes

Gel must be applied to the reusable pad before use.

Part # 8028-000014 (Case of 5)



ACCESSORIES AND SERVICES AT YOUR FINGERTIPS

Want to know more about ZOLL AED 3 accessories and services?

Contact your local authorized ZOLL distributor or visit ZOLL at
zoll.com/products/aeds/aeds-for-public-access/zoll-aed-3-for-public-access

Authorized ZOLL distributors can help you with your complete AED program, including:

- Program management and maintenance
- Compliance with local regulations and registries
- Training for your volunteer responders
- Replenishment of accessories
- Adding more AEDs to your program

ZOLL MEDICAL CORPORATION

An Asahi Kasei Group Company | 269 Mill Road | Chelmsford, MA 01824 | 978-421-9655 | 800-804-4356 | www.zoll.com

Copyright © 2020 ZOLL Medical Corporation. All rights reserved. ZOLL AED 3, Unipadz, Real CPR Help and ZOLL are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.